

HKS820

Intense Pulse Light Radio Frequency

User&Training Manual



Intense Pulse Light Safety Guidelines

The Following information is provided for the correct utilization of IPL Skin treatment system. The information includes not only the accident protection regulations the products comply with, but also the effective regarding proper use of the products.

1) Electric safety regulation

The EMC Directive 89/336/EEC amended by 92/31/EC & 93/68/EEC EN 55015:2000+A1+A2/LUM200704025/EMC EN 61000-3-2:2000+A2/LUM200704025/EMC EN 61000-3-3:1995+A1+A2/LUM200704025/EMC EN 61547:1997+A1/LUM200704025/EMC

2) Material ROHS test 2002/95/EC

Restriction of use of certain Hazardous

Substance in Electrical & electronics Equipment/ROHS.

3) Electromagnetic radiation safety regulation LVD/EN 60601-1:2006/AGC10430904GZ05E3EMC/EN 60601-1-2:2001+A1:2006/AGC10430904GZ05E1

4) ISO 13485:2003 MEDICAL/33021200005459

5) ITC MEDICAL CE

IPL SKIN TREATMENT SYSTEM belongs to II a according to MDD directive, It designed according to EN 980:2003; EN 1041:1998; EN ISO 14971:2007; EN 60601-1-2:2007; EN 61000-3-2:2006; EN 61000-3-3:1995+A1:2001+A2:2006; EN 60601-2-22:1996; EN 60601-1-4-1996+A1-1999; EN 60601-1:2006.

HKS820

Intense Pulse Light RadioFrequency

User&Training Manual

Certificate



HKS820

Intense Pulse Light
RadioFrequency
User & Training Manual

Question

Why is IPL better than other treatments?

In addition to treating multiple symptoms and locations in each session, IPL is gentle, non-invasive and safe. Because of this, there is no "down time." You see results in the complexion and the skin texture of the entire face. You can return to your normal routines and activities immediately. In fact, many people even schedule their sessions during their lunch break. Additionally, IPL treatments are customized to your specific skin type and color, providing the best possible results.

What areas can be treated with IPL systems?

The IPL System can remove hair from all parts of the body including the face, legs, arms, underarms, and back. It can also treat sensitive areas like the chest, nipples and bikini line.

What do IPL treatments feel like?

As the name implies, the IPL hand-piece cools and soothes as it glides along your skin. When the pulse of light is delivered, some patients experience a mild pinching or stinging sensation. No local anesthesia or pain medication is typically required, however some patients may prefer to use a topical anesthetic when sensitive areas are treated.

How many treatments will I need?

Most areas require 3-6 treatments, but the number of IPL treatments you may require for optimal long-term benefits depends on various factors such as your hair thickness, the area treated, and your hair growth

HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

Question

cycle. Interestingly, not all hairs are actively growing at the same time, but actually go through three distinct phases; growth, regression, and resting. IPL systems work by disabling the hair that is principally in the active growth phase at the time of treatment. Since other hairs will enter this active growth phase at different times, additional treatments may be necessary to disable all of the hair follicles in a given area.

How is IPL used to treat benign vascular lesions?
Aging, sun exposure and other factors can cause the appearance of broken capillaries and blood vessels on the face known as telangiectasias, which appear as spider veins. Sometimes broken capillaries can cause a diffused redness in the skin, called erythema, a source of social embarrassment for many people. Rosacea is a common skin condition that is characterized by diffused erythema. All of these conditions

HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

1.Product Information

1.PRODUCT INFORMATION

APPLICATION

- *Hair removal
- *Skin rejuvenation
- *red face, red nose, couperosis, spiderveins.
- *age spots, sun spots, pigmentation troubles
- *Acne
- *Skin tightening and lifting
- *Wrinkle reduction
- *Systemic collagen regeneration

FEATURES

- *1handle
- Big handle:15*50mm with 5filter
- 640-1200nm for hair removal
- 590-1200nm for red face, red nose, couperosis, spiderveins.
- 560-1200nm for Skin rejuvenation
- 430-1200nm for Acne **optional**
- 480-1200nm for age spots, sun spots, pigmentation troubles **optional**
- 750-1200nm for hair removal in Dark Skin **optional**
- *Lamp longevity:60000shots
- *Energy:15-50J
- *Pulse Number:1-15adjust
- *Pulse Width:1.0-20.0
- *Pulse Delay:1-50
- *Shot period:2-5second adjust
- *Display:8.0 inch Touch display
- *Voltage:110V or 220V
- *Power:1000W
- *Cooling system:Wind+water+semiconductor
- *RF Frequency:1/3/10MHz **Variable frequency**
- *RF Power:80W
- *RF Probe:2handle with 3probe
- *Pakcge size:66*66*110cm
- *Weight:56kg



HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

2.Installation Alignment and operation procedures


2.1 Open the Package and checking all parts.




2.2 put the handle holder on machine by screw.then connect well the handle in back of machine.



2.2 put the water to Input,when water is out of overflow means the water tank is full

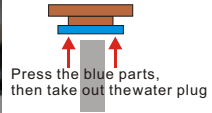

2.3 Put the filter in handle

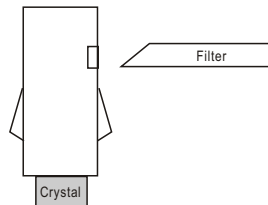


HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

3. Operate the Device

3.1 Auto interface

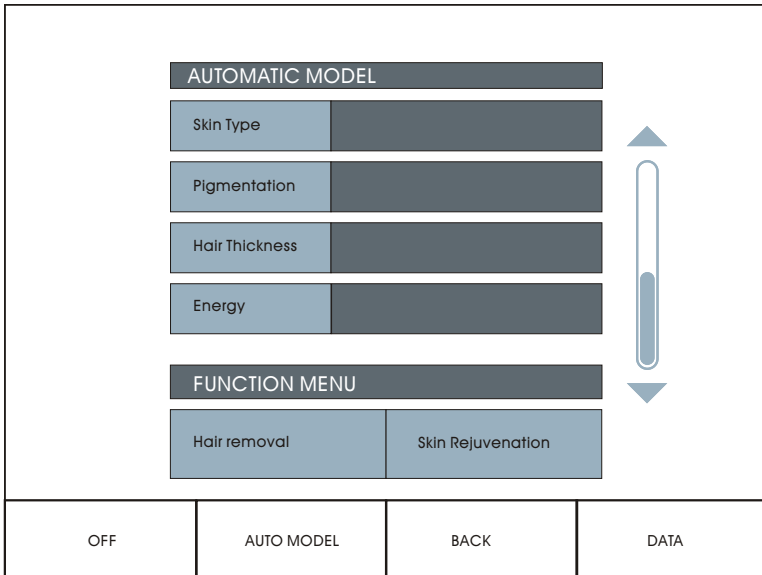

Skin Type: I -II White Skin III IV Dark Skin

Pigmentation:None Light Medium Intense

Hair Thickness:Thin Medium Thick

HKS820

Intense Pulse Light RadioFrequency

User&Training Manual

3. Operate the Device

3.2 Manual interface

<div style="border: 1px solid black; padding: 5px; margin: 10px auto; width: 80%;"> <div style="border: 1px solid black; padding: 2px; width: 60%; margin: 5px auto;">15*50mm</div> </div>	<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th colspan="2" style="text-align: left; background-color: #333; color: white;">TREATMENT</th> <th colspan="2" style="text-align: left; background-color: #333; color: white;">PARAMETER</th> </tr> </thead> <tbody> <tr> <td style="background-color: #333; color: white;">Pulse Number</td> <td style="width: 50px;"></td> <td style="background-color: #333; color: white;">15</td> <td style="width: 20px;"></td> </tr> <tr> <td style="background-color: #333; color: white;">Pulse Width</td> <td></td> <td style="background-color: #333; color: white;">20</td> <td></td> </tr> <tr> <td style="background-color: #333; color: white;">Pulse Delay</td> <td></td> <td style="background-color: #333; color: white;">50</td> <td></td> </tr> <tr> <td style="background-color: #333; color: white;">Release</td> <td></td> <td style="background-color: #333; color: white;">5</td> <td></td> </tr> <tr> <td style="background-color: #333; color: white;">Energy</td> <td></td> <td style="background-color: #333; color: white; font-size: 1.2em;">50</td> <td></td> </tr> </tbody> </table>	TREATMENT		PARAMETER		Pulse Number		15		Pulse Width		20		Pulse Delay		50		Release		5		Energy		50	
TREATMENT		PARAMETER																							
Pulse Number		15																							
Pulse Width		20																							
Pulse Delay		50																							
Release		5																							
Energy		50																							
ON	MANUAL MODEL	SAVE	DATA																						

Save: Save the parameter for current
 Release: time for each shot distance

3.3 Data interface

<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td style="background-color: #333; color: white; padding: 2px;">15*50mm</td> <td style="width: 50px;"></td> </tr> <tr> <td style="background-color: #333; color: white; padding: 2px;">Total Counter*</td> <td></td> </tr> </table>	15*50mm		Total Counter*		<table border="1" style="width: 100%; border-collapse: collapse; text-align: center;"> <tr> <td style="background-color: #333; color: white; padding: 5px;">7</td> <td style="background-color: #333; color: white; padding: 5px;">8</td> <td style="background-color: #333; color: white; padding: 5px;">9</td> </tr> <tr> <td style="background-color: #333; color: white; padding: 5px;">4</td> <td style="background-color: #333; color: white; padding: 5px;">5</td> <td style="background-color: #333; color: white; padding: 5px;">6</td> </tr> <tr> <td style="background-color: #333; color: white; padding: 5px;">1</td> <td style="background-color: #333; color: white; padding: 5px;">2</td> <td style="background-color: #333; color: white; padding: 5px;">3</td> </tr> <tr> <td style="background-color: #333; color: white; padding: 5px;">0</td> <td style="background-color: #333; color: white; padding: 5px;">Del</td> <td style="background-color: #333; color: white; padding: 5px;">Enter</td> </tr> </table>	7	8	9	4	5	6	1	2	3	0	Del	Enter
15*50mm																	
Total Counter*																	
7	8	9															
4	5	6															
1	2	3															
0	Del	Enter															
	<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td style="width: 30%;"></td> <td style="text-align: center;">▲</td> <td style="text-align: center;">▼</td> <td style="text-align: center;">BACK</td> </tr> </table>		▲	▼	BACK												
	▲	▼	BACK														

15*50mm: current shot account, press it will be 0
 total counter: all shot account, need input password for 0


HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

3. Operate the Device


3.3 Operating&setting

Before turn on machine,
must take out the water
Plug from machine back.

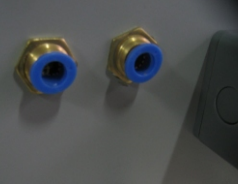

Choice one model language in Screen. The system will enter Auto model. In Auto model, the system have memory right parameter for many treatment. Only need chice right Skin type, Pigment, Hair thickness, then press the READY, the Machine can start to work. Manual model, User can input right parameter to fix sufferer treatment. Choice different filter to make different treatment. User Data show the current shot amount and total amount.

Delete total shot amount need password.

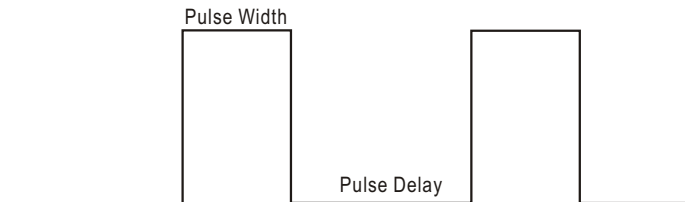

3.4 Finish the treatment.

Frist need press the ON in screen, make the lamp off. Then turn off the key switch+urgency button+main electric switch. Do not let the lamp be light when the user leave for long time.

HKS820

**Intense Pulse Light
RadioFrequency**
User&Training Manual

4.Treatment Parameter

Dark hair by lower width

Light hair by higher width

thickness hair by higher delay and lower energy

thin hair by lower delay and higher energy

dark skin by higher delay and lower energy

light skin by lower delay and high energy

When add the width that can not reduce the energy

when reduce the width that can not add the energy

SKIN TYPE	A	I - II
	B	III - IV
PIGMENTATION	1	none
	2	light
	3	medium
	4	intense
HAIR THICKNESS	C	thin
	D	medium
	E	thick

Skin rejuvenation				
	Pulse Number	Pulse Width	Pulse Delay	Energy
A1	1	9.7	0	28
A2	1	9.2	0	27
A3	1	8.8	0	26
A4	1	8.4	0	25
B1	2	9.2	10	26
B2	2	8.8	10	25
B3	2	8.4	15	24
B4	2	8.2	15	23

HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

4.Treatment Parameter

Hair Removal				
	Pulse Number	Pulse Width	Pulse Delay	Energy
A1C	1	9.5	0	36
A1D	1	8.5	0	35
A1E	1	7.5	0	33
A2C	1	8.5	0	34
A2D	1	7.5	0	33
A2E	1	6.5	0	32
A3C	2	7.5	10	32
A3D	2	6.5	10	31
A3E	2	6	10	30
A4C	3	7	15	30
A4D	3	6	15	29
A4E	3	5.5	15	28
B1C	3	9	5	35
B1D	3	8	5	34
B1E	3	7	5	33
B2C	3	8.5	10	33
B2D	3	7.5	10	32
B2E	3	6.5	10	31
B3C	3	7.5	15	31
B3D	3	6.5	15	30
B3E	3	5.5	15	29
B4C	5	6.5	17	29
B4D	5	6	17	28
B4E	5	5.5	17	27

HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

5.Treatment Procedure

5.1 Shave the hair on skin and left hair be 1-2mm.



5.2 Put the Gel on skin,Gel Thickness is around 1mm.



5.3 Touch well the crystal on skin.



5.4 After treatment, clean the Skin by warm water.

6.IPL Treatment General Knowledge

6.1.Curative effect

After first one or two treatment,the hair can lost 80%.

After 4-5 treatment,the hair will lost 95%.then no hair grow or less thin hair grow.(Acceptable)

6.2.treatment Caution And the principle

For treatment is that the patient should feel little Acceptable pain.Keep the skin be cool.

Never shot to eyes

6.3. Taboos

- a.To avoid using the machine on the abdomen of the pregnant women or women in menses
- b.To avoid using the machine on around the parts under which metal,plastic or silicon we buried.
- c.To avoid using the machine on around the medical and electric instrument,such as the hearing aid,heart starter, mechanical heart,etc.
- d.To avoid using the machine on the epileptic and those in bad health
- e.To avoid using the machine on the wound around the treatment area and the areas with metals inside of the body.
- f.To avoid using the machine on people who has high blood pressure or diabetes or serious skin diseases or heart diseases, or who is sensitive to electric current,or who is receiving hormone treatment.

Radio Frequency Safety Guidelines

The Following information is provided for the correct utilization of IPL Skin treatment system.the Infirmination includes not only the accident protection regulations the products comply with, but also the effective prearding proper use of the products.

1)Electric safety regulation

The EMC Directive 89/336/EEC amended by 92/31/EC & 93/68/EEC EN 55015:2000+A1+A2/LUM200704025/ EMCEN 61000-3-2:2000+A2/LUM200704025/ EMCEN 61000-3-3:1995+A1+A2/LUM200704025/ EMCEN 61547:1997+A1/LUM200704025/EMC

2)Material ROHS test 2002/95/EC

Restriction of use of certain Hazardous

Substance in Electrical&electronics Equipment/ROHS.

3)Electromagnetic radiation safety regulation LVD/EN60601-

1:2006/AGC10430904GZ05E3EMC/EN60601-1-2:2001+A1:2006/AGC10430904GZ05E1

4)ISO13485:2003 MEDICAL

HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

Radio Frequency Safety Guidelines

What is Radio Frequency?

RF is used for dermal skin heating with sparing of the epidermis. This will result in non-ablative skin tightening and lifting without injection and surgery.

What is the advantage?

The results are a firmer, refreshed and more youthful face, lifted brows and eyelids, less wrinkles and sharper jaw line again. Eyebags and double chin will improve with this treatment too.

How does it work?

First, heat generated by RF waves shrinks the bundle of collagen under the skin. The collagen becomes shorter and the skin tightens. The second phase begins as the body responds. The collagen will be gradually healed, renewed and even more collagen will be produced. The skin tightens and overall look is improved.

What do I expect during the treatment?

gentle warming sensation is felt throughout the procedures. As there is minimal redness after the treatment, normal activity can be resumed immediately.

How many sessions are required?

1times/2weeks

5-6 sessions is recommended for optimal longer lasting result in treating the wrinkles and lifting the face.




HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

7.RF Operate the Device

7.1 Operate interface

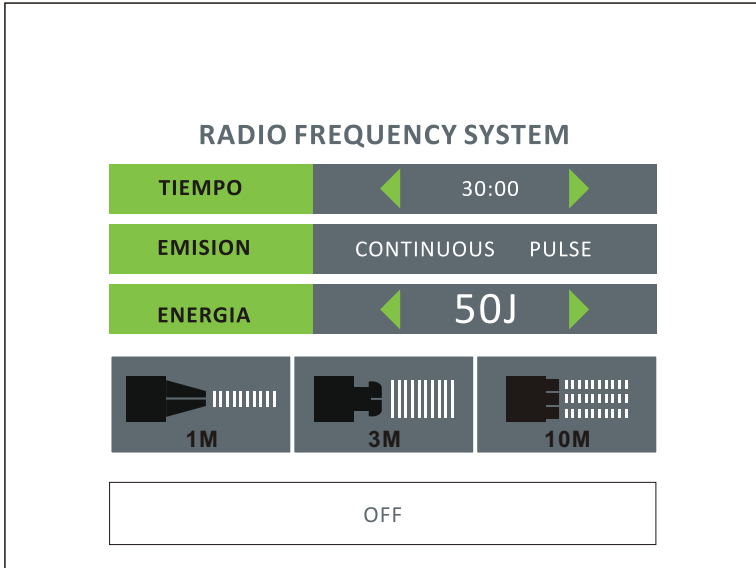

HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

7.RF Operate the Device


7.1 Operate interface



Small Probe forFace

Big Probe forbody/neck

Tripolar Forbody/Jowl

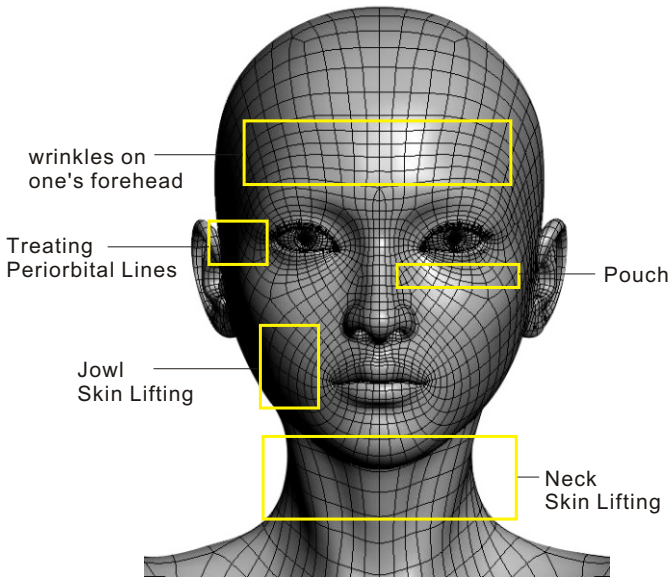








HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

7.RF Operate the Device

7.2 Manual Model Treatment

Choice right Probe as you want and adjust the timer and energy,press OFF be ON,the Device start to work.
User need press the switch in handle.

	Adjust Button
	Timer Adjust,Max is 30Minutes
	Continuous or discontinuous
	Energy Adjust
	Choice small bipolar probe 1MHz
	Choice big bipolar probe 3MHz
	Choice tripolar probe 10MHz


HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual


8.Treatment Procedure

8.1 .Operateion technigue:Let the treatment probe touch the target area ad mover(Or Slide)from the inner to the outer and round and round and from the lower place to the higher, when you do the treatment ont he upper lid,Pls move a bit more quickly,when the gel on the face get drier,please add some GEL and go on with the treatment.

8.2.Face(Includeing cheek,eye,forhead,lip,neck)
Operation procedure(About 50minutes for whole face,15-20minutes for neck)

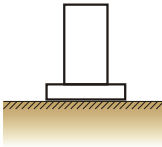

- a. carry out a skin evaluation before treatment.
- B. Clean skin
- c. tone skin
- d. massage with GEL
- e. clean skin
- f. smear gel on skin and apply treatment with trobe.
- g. let the probe touch the skin and make treatment

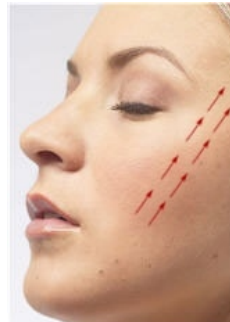

Course of treatment:One course of treatments includes 5treatments at an interval of every 10-15days.

Curative effect

After one treatment,the heated subcutaneous fiber shrinks instantly,therefore fins wrinkles are improved,and skin becomes tender,sommth and compact.After onecourse of treatment:the collagen protein is regenerated for long term,skin becomes tender and smooth,wrikles are reduced, and skin is tightened and compact.

HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

8.Treatment Procedure

8.3.Body(striae gravidarum)

Operation procedure(About 60-90minutes for topical area)

a.Carry out a skin evaluation before treatment.

B.clean skin

c.massage with body slimming distillate

d.Smear body slimming diatillate on skin and apply treatment with probe.


e.clean skin

f.let the probe touch the skin and make treatment.

Course of treatment:One course of treatments includes 5 treatments at a interval of every 10-15days.

Curative effectAfter one treatment,the heated subcutaneous fiber shrinks instantly,therefore skin become smooth,compact and elastic.

After one course of treatment;Skin becomes compact and flabby skin is improved.

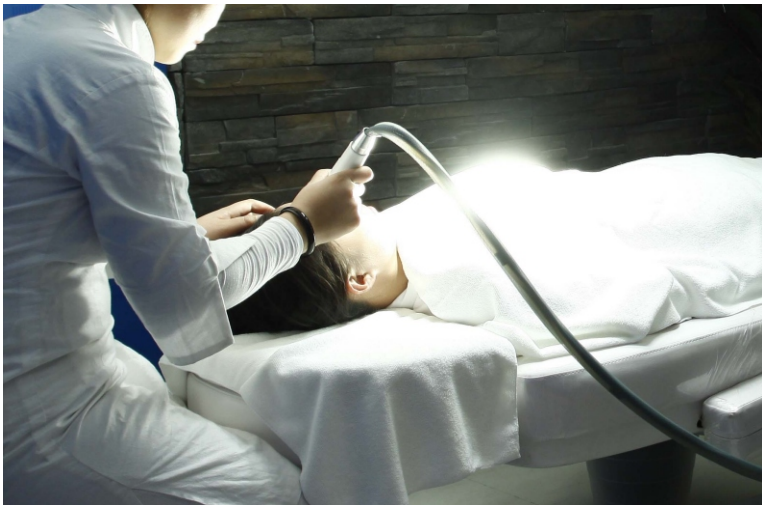

HKS820

Intense Pulse Light
RadioFrequency
User&Training Manual

9.Treatment Caution

9.1.for treatment is that the patient should feel needling&hot.

9.2.The treatment head must be kept on the skin (sliding); otherwise, the skin will get burned.Before the treatment head leaves the skin, you must loose the button on the head or loose the foot switch.

9.3.When the temperature of the treatment head is over too high, you should stop treatment.The energy you adopt should not be too high; otherwise, the skin will get burned.

9.4.Please carry out skin moisturizing and sun proof cream work after treatments, do not wash faces with hot water (use water under the body temperature), and do not go for hot spring and massage.

10.Treatment Taboos

10.1.The abdomen of the pregnant women or women in menses

10.2.To avoid using the machine on around the parts under which metal,plastic or silicon wee buried.

10.3.To avoid using the machine on the epileptic and those in bad health

10.4.To avoid using the machine on people who has high blood pressure or diabetes or serious skin diseases or heart diseases, or who is sensitive to electric current,or who is reveiving hormone treatment.

HKS820

Intense Pulse Light RadioFrequency

User&Training Manual

11. Treatment Parameter

Therapy Position	Probe	Energy	Mode	Timer
eye side	small	15-20	■	10
	big	22-28	■	5
forehead	small	14-17	■	8
	big	22-28	■ ■ ■ ■	5
jowl	Tripolar	20-26	■	10
neck	big	22-28	■	15
Abdomen	big	22-28	■	10
buttocks	Tripolar	22-28	■ ■ ■ ■	10
Low-necked	Tripolar	18-24	■	10
arm	Tripolar	22-28	■ ■ ■ ■	10
leg	Tripolar	22-28	■ ■ ■ ■	10

12. HKS817C RF Device Superiority

- 12.1. Whole ABS Plastic Mould.
- 12.2. Variable frequencies depending on the chosen Treatment.
- 12.3. Frequency of steady emission.
- 12.4. It does not present any contagious risk.
- 12.5. Maximum security for the client, default programs
- 12.6. Permits to adapt the used power used according to the zone and deepness to work.



We Are Looking for Sale agent over the world.